


# UniTip APPLICATORS

POLYSACCHARIDE HEMOSTATIC POWDER  
APPLICATION SYSTEM FAMILY

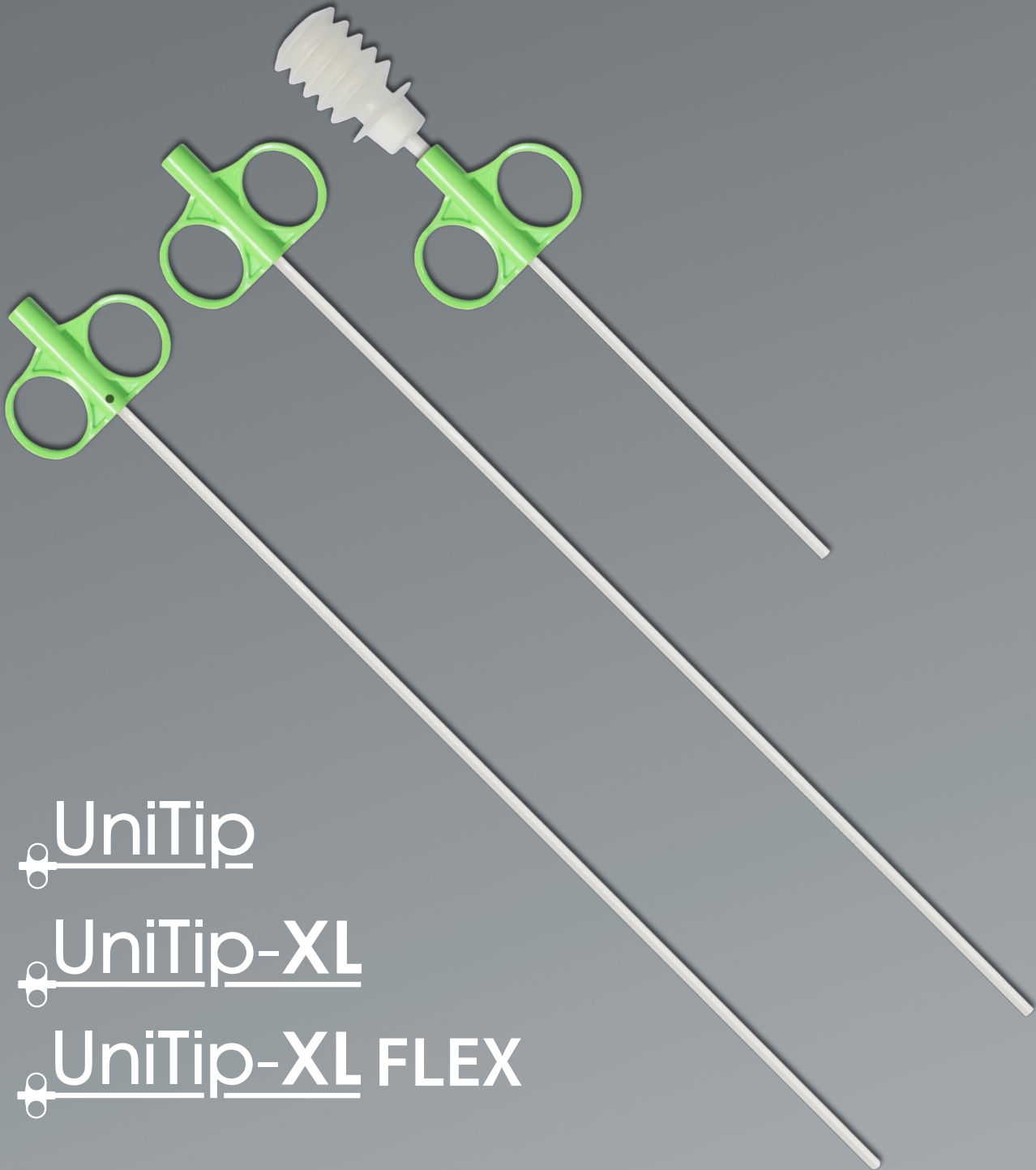

UniTip

UniTip-XL

UniTip-XL FLEX

**3**

**UNITIP APPLICATORS DESIGNED FOR  
USE WITH STARSIL HEMOSTAT®**

FOR LAPAROSCOPIC, MINIMALLY INVASIVE, ENT & GENERAL SURGERY

# THE UNITIP APPLICATOR FAMILY IS DESIGNED FOR THE EASY AND FAST APPLICATION OF STARSIL HEMOSTAT® HEMOSTATIC POWDER

STARSIL HEMOSTAT® together with the UniTip applicators is suitable for almost all kinds of interventions, surgical procedures and wounds. More and more surgeries and interventions are being carried out minimally invasively. The UniTip applicators now enable STARSIL HEMOSTAT® to be used even more flexibly.


The easy plug-in mechanism guarantees a secure airtight connection. The ergonomically shaped handpieces of the applicators ensure comfortable use, and the atraumatic applicator tips protect the surrounding tissue from injuries. Significantly improving the handling and enabling safe application of the STARSIL HEMOSTAT® powder even in difficult-to-reach regions.



## INNOVATIVE NANO STRUCTURED MATERIAL WITH LOTUS EFFECT

The innovative nanostructured surface of the applicator tubes with the lotus effect allows the application of the STARSIL HEMOSTAT® topical hemostat without the need for a hand pump or external compressed air source and minimizes the sticking of blood or other liquids.

The UniTip Applicator Family offers three options, providing the right length and diameter for Laparoscopic, Robot-Assisted, ENT, and General surgeries. The opening of the UniTip (14 cm) and UniTip XL applicators is located at the front of the applicator tip. With the new UniTip XL Flex (38 cm) applicator, the aperture is positioned on the side of the applicator tip, facilitating finer control over the dispersion of the STARSIL HEMOSTAT® powder, particularly in scenarios necessitating enhanced guidance of the spray direction.

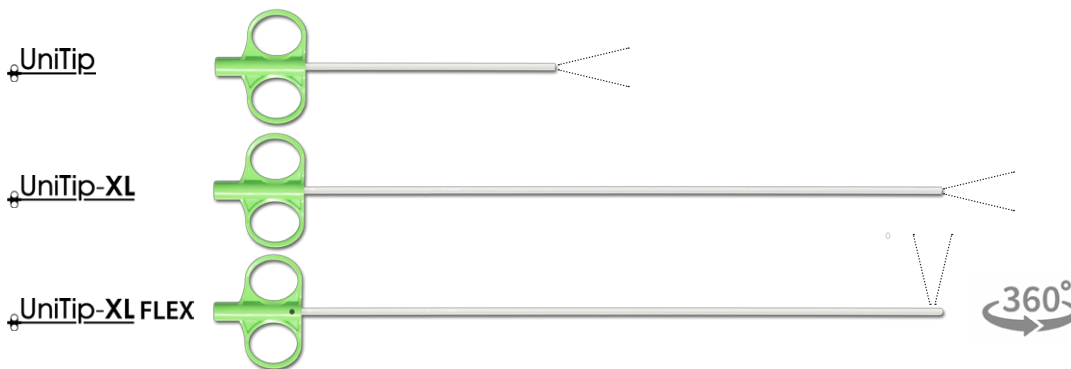

Order no.	Product	Diameter	Length	Sales unit
UT0014	UniTip Applicator	3,0 mm	14 cm	10 pcs.
UT0038	UniTip XL Laparoscopic Applicator	3,0 mm	38 cm	10 pcs.
UTF038	UniTip XL Flex Laparoscopic Applicator	3,0 mm	38 cm	10 pcs.

DISTRIBUTED BY:

**HEMOSTAT  
MEDICAL**  
HEMOSTAT MEDICAL GMBH  
BECKELMANNSWEG 10 · D-46342 VELEN · GERMANY  
TEL +49 (0) 2863 9246446 · FAX +49 (0) 2863 2450  
info@hemostatmedical.com · www.hemostatmedical.com